

HYPER-ANGULATED VL




Blade

Size 3 (preferred)
Lubricate inner curvature



Stylet

Straight to cuff then 70°
Do NOT lubricate
Reverse loaded into ETT

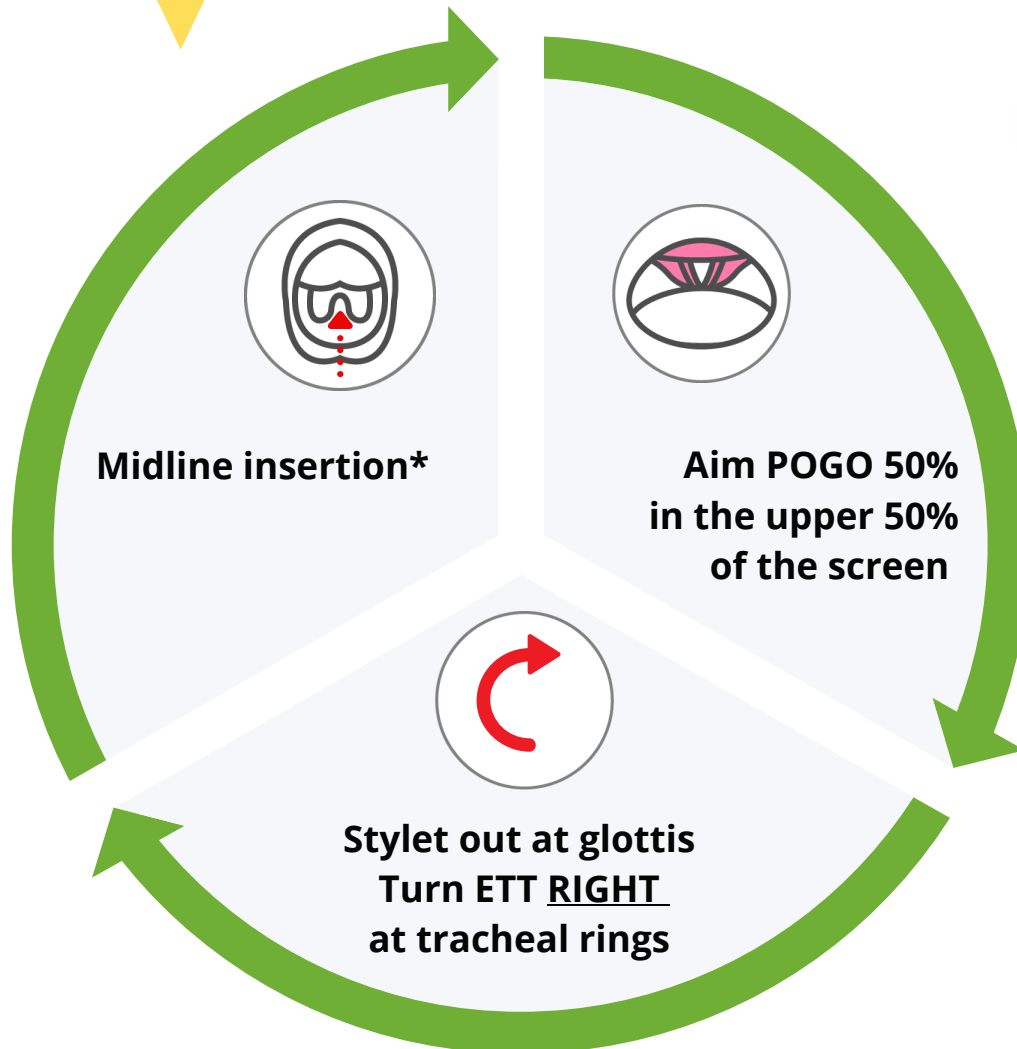

***If difficulty inserting the blade**

Insert blade at 90 degrees rotation then re-orientate midline

Difficult ETT delivery to cords

- Lift laryngoscope + stylet in direction of ceiling.
- Neck flexion
- Bougie (e.g. Flexi-tip)

Difficult advancing ETT

- Advance ETT after "popping off" stylet by 3-5cm
- Reinforced or Fast-trach ETT